



<https://happysield.github.io>

Headband

Elastic Strap
(20-30x350)mm

Strap Holder

Transparent Plastic sheet
(30x33x0.5) mm

Shield

Transparent Plastic sheet
(297x315x0.5) mm

Second fold


First fold

Flat







<https://happysield.github.io>

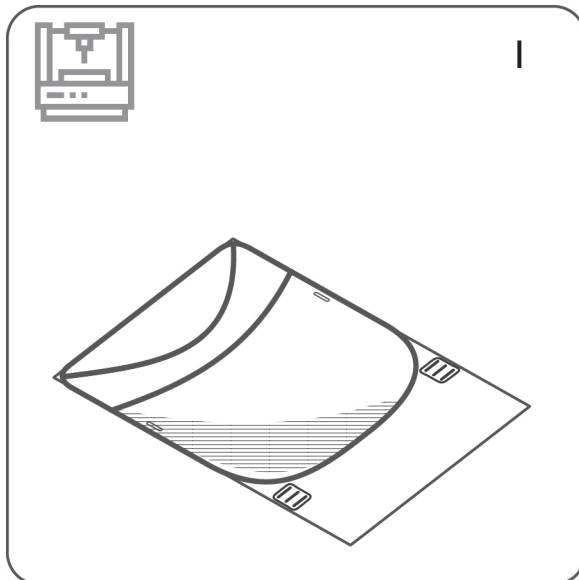


Download and use our laser cutting template. Place the plastic sheet and laser cut. The mask cutting area is 297x320mm. You can use A3 sheets if large plastic sheets are not available.

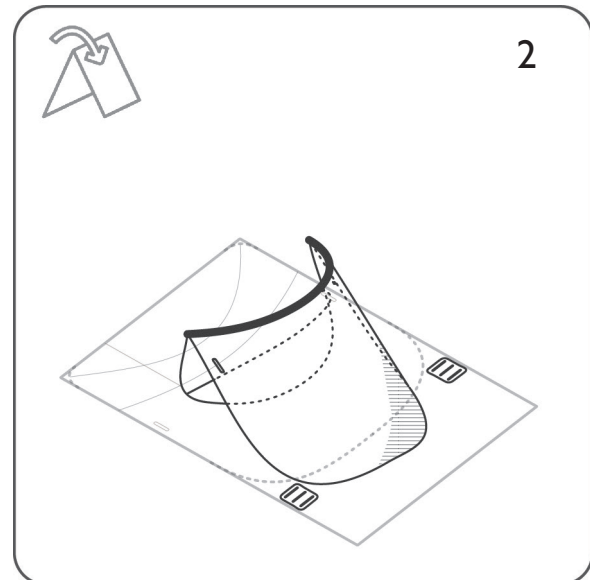

To fold the first curve, start with one corner; if some segments do not fold, score again with a pen and ruler.

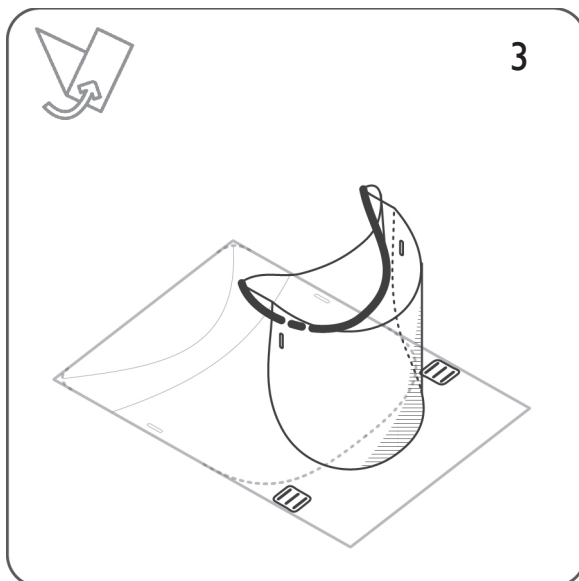

Fold the second curve.

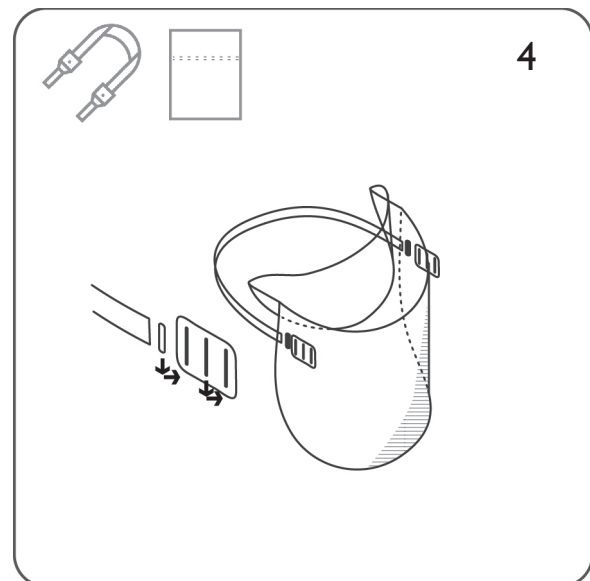

Fix the strap using the holders you cut from the template. Clean and use or clean and place in a ziploc bag for shipping.